



# LX-Mask Plus

Red & Blue Light Therapy Face Mask

# User Manual

## 1. About American Wellness Authority:

Our journey began long ago, first at a retail store in North Carolina and then online. Over the years, we've sold some of the country's most iconic brands. Ranging in categories like fitness accessories, professional beauty, healthcare, gourmet food, personal care, and many in between. We've sold thousands of SKUs to happy customers all over North America, and we feel truly blessed for the opportunity.


During the Pandemic, we realized how difficult it was to continue delivering our best sellers. We couldn't replenish our stock and there wasn't much we could do. We didn't have control of the production. We couldn't change any designs or product capabilities. With our options limited, we decided to double down on our knowledge and experience, and bring to you American Wellness Authority, or AWA for short.

Since then, our team of 10 spent countless months researching and debating. Searching for ways to deliver practical products, one's we'd be proud of. Our focus was on things that make people feel good and could enhance their lives in one way or another.

So we set out on our journey, and here we are today. We couldn't be more excited to share this with you. And we sincerely appreciate your taking an interest in our product and our company. We're a small, family-run business and there isn't enough we can do to say, thank you.

## 2. About Light Therapy:

Light is a composition of various electromagnetic waves. The white light that we see has a spectrum of 7 different color lights with different wavelengths.

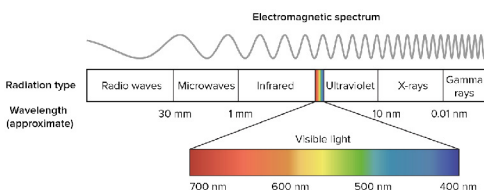

Different wavelengths of visible light penetrate the skin at different depths. Depending on how deeply they penetrate, LED lights to have different biological effects that treat various skin and muscle-related

### Red Light Therapy (630nm):

It can improve cell activity, speed up metabolism, and promote collagen formation. It has the effects of whitening and lightening, rejuvenating & wrinkling, increasing skin elasticity, improving skin condition, repairing damaged skin, reducing pores and tightening skin, anti-aging and anti-oxidation.

### Blue Light Therapy (470nm):

It has the effect of bactericidal, anti-inflammatory and inhibiting inflammation, and can remove acne and sores, sterilize and reduce inflammation, repair skin, and leave no scars and pigments.

### Green Light Therapy (520nm):

It can reduce the secretion of skin oil, balance the skin, effectively relieve stress, dredge lymph and drain swelling.

### Purple Light Therapy (390nm):

It is a dual-frequency light of red light and blue light. It combines the effects of two kinds of phototherapy, especially in the treatment of acne and acne scars and has a particularly good effect and repairing effect.

### Cyan Light Therapy (490nm):

It can increase cell energy step by step, and has a good effect on metabolism.

### Yellow Light Therapy (590nm):

Replenish energy to skin cells, decompose pigments, promote lymphatic drainage, improve skin roughness and wrinkles, treat skin redness, fever, ringworm, etc., and enhance immunity.

### White Light Therapy (mix):

Penetrates the deep layers of the skin, accelerates skin metabolism, breaks down pigment spots, improves fine lines and sagging skin.

## 3. Direction of Use:

- Preparation:** Clean and dry your face. Avoid using coarse towels that may lead to rough skin.
- Set-Up:** Connect the mask and power control box. Ensure a secure connection.
- Power On:** Press the "Color Select/Start" button to turn on the device. Default time is 30 minutes.
- Mode Selection:** Use the "Color Select/Start" button to choose the desired light mode.
- Adjustment:** Use "Time -" and "Time +" buttons to regulate working time (0-60 minutes) as per your preference.
- Caution:** Avoid prolonged staring at the light source for eye safety.
- Brightness and Flicker:** Adjust light intensity (01 - 09) and flashing mode using respective buttons.
- Micro Electricity:** Customize micro electric strength (01 - 09) as needed.
- Wearing the Mask:** Position the mask comfortably on your face, adjusting the magic band for a secure fit.
- Using Skincare Products:** Apply skincare products evenly before using the mask for enhanced benefits.
- Completing the Session:** After the set time, the power control box will automatically turn off the LED photons. Remove the mask and proceed with cleaning.

## 4. Function & Modes:

The LX-Face Mask Plus is designed to provide the following functions:

- Increasing Cell Activity:** Promoting cell metabolism for enhanced skin health.
- Stimulating Collagen and Fibrous Tissue Production:** Filling the skin for improved elasticity.
- Accelerating Blood Circulation:** Enhancing skin tone and reducing dullness.
- Micro Electricity Massager:** Provides gentle facial massages to rejuvenate and refresh the skin.



### Control The Device Remotely



## 5. Technical Specifications:

Material	ABS+LED Lights	Color	White
Actual Power	36W	Wavelength	
LED Quantities	192	Product Size	12 × 12.2 × 9.1 in
Output Current		Input Voltage	
Single pcs Weight		Temperature	
Package Dimension	16.9 × 11.8 × 3.1 in	Product Weight	3.3 lbs
Light Colors	Red, Purple, Green, Yellow, Blue, White, Cyan		

## 6. Precautions:

---

- **Not Waterproof:** The mask is not waterproof. Avoid contact with water to prevent damage.
- **Not Suitable for Children and Pregnant Women:** This product is not intended for use by children or pregnant women.
- **Use Maximum 30 Minutes:** Limit usage to a maximum of 30 minutes per session.
- **Do Not Stare at Light:** Avoid direct eye contact with the LED lights to protect your eyes.
- **Not for Use While Sleeping:** Do not use the mask while sleeping or in a reclined position.
- **Keep Away from Children:** Store the mask in a safe place out of reach of children.

## 7. Cleaning:

---

- Gently wipe the mask with a soft, damp cloth after each use. Avoid using harsh chemicals.
- Ensure the micro electricity massager tips are cleaned thoroughly.
- Do not immerse the mask in water or use harsh chemicals for cleaning.

## 8. Storage:

---

- Store the mask and any additional components in a cool, dry place away from direct sunlight.

## 9. FAQs:

---

### What are the EMF levels?

For best results, use the mask 3-4 times per week.

### Can I use skincare products with the mask?

Yes, you can apply your usual skincare products before using the mask.

### Can I share the mask with others?

It is recommended to use the mask individually for hygiene purposes.

### For How long should the device be used for effective treatments?

Commit to a three-week usage plan. While many people achieve relief from Red Light Therapy rather quickly, others see measurable relief incrementally over time.

Therapy response time varies depending on skin type.

### Does light therapy work the same for all skin types?

Light therapy is an effective option for all skin types and pigmentation. It's just important to use the treatment regularly to achieve optimal results.

### Does red light therapy work through clothes, bandages, or dressings?

Unfortunately, red light therapy does not work through clothes, bandages, or dressings. This is because the red light must glow on your skin for it to be absorbed correctly. If used on clothes, the light will be absorbed by the fabric.

## 10. Warranty:

---

For additional information, feedback, and warranty registration, please don't hesitate to contact us. Our dedicated support team is here to assist you.

### Contact Information

*Email:* [Contact@awastore.us](mailto:Contact@awastore.us)

### Warranty Details:

Enjoy peace of mind with a 1-year warranty for your LX-Face Mask -Red & Blue Light Therapy Face Mask

### Warranty Registration:

To register your product for warranty coverage and access exclusive benefits, simply scan the QR code below or visit our website.



We value your satisfaction and are committed to providing

## 11. Disclaimer:

---

This device has not been evaluated by the U.S. Food & Drug Administration (FDA). Our products are designed to maintain or encourage a healthy lifestyle and are not intended to diagnose, treat, cure, or prevent any disease. Our products are low-risk, general wellness/fitness products that do not require FDA clearance, in accordance with the "General Wellness: Policy on Low-Risk Devices" draft released on January 20, 2015.

## 12. Contact Us:


---

You can contact the AWA team via Email: [contact@awastore.us](mailto:contact@awastore.us)  
You can visit [www.awastore.us](http://www.awastore.us) for more information.


 American Wellness Authority

 American Wellness Authority

 [awa\\_products](#)

 [awastore.us](#)

 [awastore\\_us](#)

 [awastore\\_us](#)