



LX-Cap+
Red Light Therapy Cap
User Manual

1. About American Wellness Authority:

Our journey began long ago, first at a retail store in North Carolina and then online. Over the years, we've sold some of the country's most iconic brands. Ranging in categories like fitness accessories, professional beauty, healthcare, gourmet food, personal care, and many in between. We've sold thousands of SKUs to happy customers all over North America, and we feel truly blessed for the opportunity.


During the Pandemic, we realized how difficult it was to continue delivering our best sellers. We couldn't replenish our stock and there wasn't much we could do. We didn't have control of the production. We couldn't change any designs or product capabilities. With our options limited, we decided to double down on our knowledge and experience, and bring to you American Wellness Authority, or AWA for short.

Since then, our team of 10 spent countless months researching and debating. Searching for ways to deliver practical products, one's we'd be proud of. Our focus was on things that make people feel good and could enhance their lives in one way or another.

So we set out on our journey, and here we are today. We couldn't be more excited to share this with you. And we sincerely appreciate your taking an interest in our product and our company. We're a small, family-run business and there isn't enough we can do to say, thank you.


2. About Light Therapy:

Light is a composition of various electromagnetic waves. The white light that we see has a spectrum of 7 different color lights with different wavelengths.

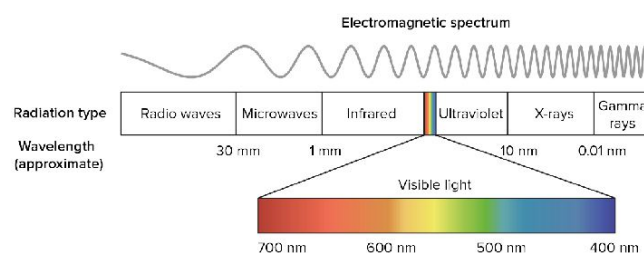


Different wavelengths of visible light penetrate the skin at different depths. Depending on how deeply they penetrate, LED lights to have different biological effects that treat various skin and muscle-related issues.

Red Light Therapy: (660 nm)

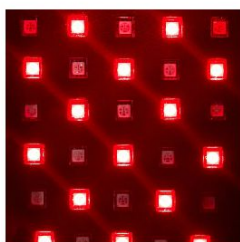

- It targets skin cells (fibroblasts), which play a role in collagen production. This protein makes up a large part of connective tissue and helps the skin to recover when it's harmed.
- Reduce skin cell inflammation and redness and improve blood circulation.
- Reduce wrinkles, and dull pigmentation, relieve skin fatigue, restore skin elasticity and luster, and promote skin rejuvenation and deep repair.

Infrared or Near-Infrared Light Therapy: (850 nm)

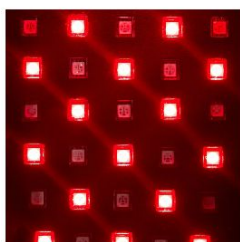

- Penetrate beyond Dermis and reach muscle tissues.
- Get absorbed by the Mitochondria (also called the Powerhouse of the cell) that accelerate the healing process in the tissues and promote blood circulation that reduces pain effectively.
- Allow skin products to absorb better by enhancing skin elasticity and moisture retention. Helps in shrinking pores.

3. Direction of Use:

1. Place the LX Cap+ on your head, ensuring it covers the entire scalp.
2. Connect the remote control to the cap.
3. Connect the remote to a power adapter in a socket or use a power bank.
4. Turn on the device using the remote control.
5. Select the desired level (low, medium, or high) using the remote.
6. Sit back, relax, and allow the LX Cap+ to work its magic for 10-30 minutes.



4. Function & Modes:

The LX Cap+ is a cutting-edge device designed to promote hair growth and scalp rejuvenation.

Its key components include:

1. Professional-grade LEDs
2. Built-in timer for convenience
3. Remote control for easy operation



5. Technical Specifications:

Material	PU+ Fabric	Shell Color	Black
Actual Power	75 W	Wavelength	850nm : 940nm
LED Quantities	150 pcs	Product Size	14.2 × 6.3 × 9.4 in
Output Current		Input Voltage	
Single pcs Weight		Temperature	
Innerbox Dimension		G.W.	
Carton Dimensions		Product Weight	0.4 lbs

6. Precautions (DOs and DON'Ts):

- **Do's**

1. Consult a doctor if you feel unwell during or after use.
2. Use this product under the guidance of a doctor if you have sensitive skin or allergic dermatitis.
3. Unplug the product's connecting wire by holding the plug and pulling it out from the socket.
4. Store the product in a ventilated and dry place after use.
5. Clean the product by wiping it gently with a soft cloth on a regular basis.

- **Don'ts**

1. Do not use the product in places with high humidity and excessively high temperature, such as bathrooms.
2. Do not disassemble, modify, immerse in water or throw the product into fire.
3. Pregnant women, cancer patients, or those who have recently undergone surgery should not use this product.
4. Do not apply heavy pressure to the product.
5. Do not soak or rinse in water, and do not use strong cleaning agents to clean.

7. Cleaning:

- Ensure the device is turned off and disconnected from any power source before cleaning.
- Take a soft, damp (not wet) cloth. You may dampen it with warm water if necessary. Ensure the cloth is not dripping with water to prevent moisture from seeping into the device.
- Carefully wipe the fabric surface of the cap with the damp cloth. Gently rub any areas that need cleaning.
- Do not soak the cap or use excessive water as it may damage the LEDs and the device's electronic components.
- After cleaning, allow the cap to air dry completely before using it again. Ensure it's completely dry to prevent any potential moisture-related issues

8. Storage:

- Store the cap in a cool, dry place when not in use.
- Ensure that it is kept away from direct sunlight, extreme temperatures, and moisture.
- Place the cap in a protective bag or cover, if available, to prevent dust accumulation.
- Keep the remote control and power cord with the cap for convenient access when needed.

9. FAQs:

Can I use the LX Cap+ on wet hair?

It is recommended not to use the cap on wet hair immediately after showering, swimming, or bathing. Ensure your hair is dry for optimal results and safety.

Can multiple family members use the same cap?

While the cap is designed for individual use, it can be sanitized between uses for hygiene purposes. However, it's recommended to have a separate cap for each user for a more personalized experience.

How do I clean and maintain the LX Cap+?

Use a soft, damp cloth to gently clean the caps's surface. Avoid using chemicals or excessive moisture. Store it in a cool, dry place when not in use.

Can the LX Cap+ be used with other hair treatments or medications?

Consult with a healthcare professional before combining the use of the cap with other treatments to ensure safety and effectiveness.

Is the LX Cap+ safe for children and pets?

Keep the cap out of reach of children and pets to prevent accidental use.

How does the LX Cap+ work, and what is its primary function?

The LX Cap+ employs red light therapy to stimulate hair follicles and promote hair growth. Its primary function is to address hair loss, thinning hair, and scalp health.

10. Warranty:

For additional information, feedback, and warranty registration, please don't hesitate to contact us. Our dedicated support team is here to assist you.

Contact Information

Email: Contact@awastore.us

Warranty Details:

Enjoy peace of mind with a 1-year warranty for your LX - Cap+ Red Light Therapy

Warranty Registration:

To register your product for warranty coverage and access exclusive benefits, simply scan the QR code below or visit our website.



We value your satisfaction and are committed to providing exceptional customer support. Your feedback and warranty registration help us ensure your product's continued performance.

11. Disclaimer:

This device has not been evaluated by the U.S. Food & Drug Administration (FDA). Our products are designed to maintain or encourage a healthy lifestyle and are not intended to diagnose, treat, cure, or prevent any disease. Our products are low-risk, general wellness/fitness products that do not require FDA clearance, in accordance with the "General Wellness: Policy on Low-Risk Devices" draft released on January 20, 2015.

12. Contact Us:


You can contact the AWA team via Email: contact@awastore.us

You can visit www.awastore.us for more information.


 American Wellness Authority

 American Wellness Authority

 [awa_products](#)

 [awastore.us](#)

 [awastore_us](#)

 [awastore_us](#)