



LX-Slippers

Red Light Therapy Panel

User Manual



1. About American Wellness Authority:

Our journey began long ago, first at a retail store in North Carolina and then online. Over the years, we've sold some of the country's most iconic brands. Ranging in categories like fitness accessories, professional beauty, healthcare, gourmet food, personal care, and many in between. We've sold thousands of SKUs to happy customers all over North America, and we feel truly blessed for the opportunity.


During the Pandemic, we realized how difficult it was to continue delivering our best sellers. We couldn't replenish our stock and there wasn't much we could do. We didn't have control of the production. We couldn't change any designs or product capabilities. With our options limited, we decided to double down on our knowledge and experience, and bring to you American Wellness Authority, or AWA for short.

Since then, our team of 10 spent countless months researching and debating. Searching for ways to deliver practical products, one's we'd be proud of. Our focus was on things that make people feel good and could enhance their lives in one way or another.

So we set out on our journey, and here we are today. We couldn't be more excited to share this with you. And we sincerely appreciate your taking an interest in our product and our company. We're a small, family-run business and there isn't enough we can do to say, thank you.


2. About Light Therapy:

Light is a composition of various electromagnetic waves. The white light that we see has a spectrum of 7 different color lights with different wavelengths.

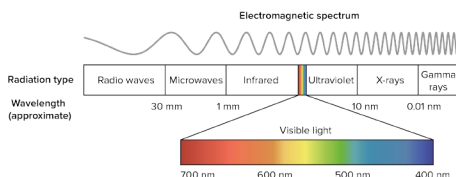


Different wavelengths of visible light penetrate the skin at different depths. Depending on how deeply they penetrate, LED lights to have different biological effects that treat various skin and muscle-related issues.

Red Light Therapy: (660 nm)

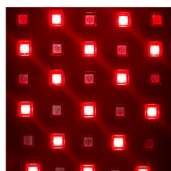

- It targets skin cells (fibroblasts), which play a role in collagen production. This protein makes up a large part of connective tissue and helps the skin to recover when it's harmed.
- Reduce skin cell inflammation and redness and improve blood circulation.
- Reduce wrinkles, and dull pigmentation, relieve skin fatigue, restore skin elasticity and luster, and promote skin rejuvenation and deep repair.

Infrared or Near-Infrared Light Therapy: (850 nm)

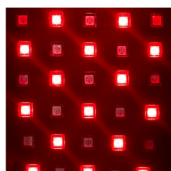


- Penetrate beyond Dermis and reach muscle tissues.
- Get absorbed by the Mitochondria (also called the Powerhouse of the cell) that accelerate the healing process in the tissues and promote blood circulation that reduces pain effectively.
- Allow skin products to absorb better by enhancing skin elasticity and moisture retention. Helps in shrinking pores.

3. Direction of Use:

Option 1: Plug in with an adapter

Option 2: Use a 5V to 12V USB cable to plug in the power bank

1. Wrap your body, choose a part of the area where you will place it, connect the adapter and the controller into the power connector.
2. Plug into an AC outlet and long press the power button to turn on the device or use a power bank to power the device.
3. According to the sensitivity of the human body to temperature, choose different levels. It is recommended to try from the lowest gear for the first time.
4. Press the appropriate button to set the time and power usage.

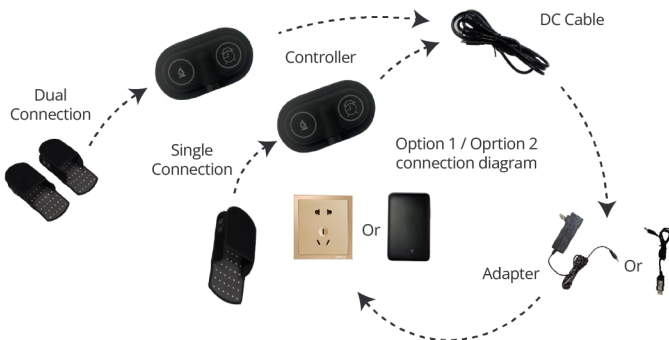


4. Function & Modes:

Controller details:

- a. Voltage input DC 12V
- b. Time gear setting: 3 gears (green, blue, red), green: 10 minutes, blue: 20 minutes, red: 30 minutes
- c. Power gear setting: 3 gears (green, blue, red), green: minimum, blue, red: maximum Controller

Operating Instructions:

- a. After the power is connected, press and hold the time setting key for 3 seconds to turn on/off
- b. Power setting key: green light (minimum), blue light, red light (maximum)
- c. Time setting key: green light 10 minutes, blue light 20 minutes, red light 30 minutes

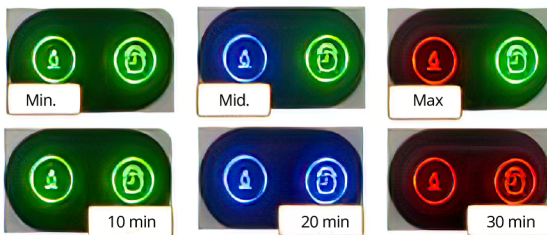

5. Technical Specifications:

Material	PU + Fabric	Color	Black
Actual Power		Wavelength	660nm:850nm
LED Quantities	90 Pcs	Product Size	12.6 × 5.1 × 2.0 in
Output Current		Input Voltage	AC100-240V
Single pcs Weight		Temperature	
Innerbox Dimension		C.W.	
Carton Dimensions	11.8 × 9.4 × 3.5 in	Product Weight	3.3 lbs

6. Precautions (DOs and DON'Ts):

● Do's

1. Consult a doctor if you feel unwell during or after use.
2. Use this product under the guidance of a doctor if you have sensitive skin or allergic dermatitis.
3. Unplug the product's connecting wire by holding the plug and pulling it out from the socket.
4. Store the product in a ventilated and dry place after use.
5. Clean the product by wiping it gently with a soft cloth on a regular basis.

● Don'ts

1. Do not use the product in places with high humidity and excessively high temperature, such as bathrooms.
2. Do not disassemble, modify, immerse in water or throw the product into fire.
3. Pregnant women, cancer patients, or those who have recently undergone surgery should not use this product.
4. Do not apply heavy pressure to the product.
5. Do not soak or rinse in water, and do not use strong cleaning agents to clean.



7. Cleaning:

- Ensure the device is turned off and disconnected from any power source before cleaning.
- Take a soft, damp (not wet) cloth. You may dampen it with warm water if necessary. Ensure the cloth is not dripping with water to prevent moisture from seeping into the device.
- Carefully wipe the fabric surface of the slippers with the damp cloth. Gently rub any areas that need cleaning.
- Do not soak the slippers or use excessive water as it may damage the LEDs and the device's electronic components.
- After cleaning, allow the slippers to air dry completely before using it again. Ensure it's completely dry to prevent any potential moisture-related issues

8. Storage:

- Store the slippers in a cool, dry place when not in use.
- Ensure that it is kept away from direct sunlight, extreme temperatures, and moisture.
- Place the slippers in a protective bag or cover, if available, to prevent dust accumulation.
- Keep the remote control and power cord with the slippers for convenient access when needed.

9. FAQs:

What is the purpose of LX-Slippers?

LX-Slippers combine the comfort of footwear with the benefits of red light therapy. They help soothe tired feet, improve blood circulation, and provide relaxation and pain relief.

How does red light therapy work in LX-Slippers?

LX-Slippers emit gentle yet effective red light, which is known to have therapeutic benefits. This light therapy helps in rejuvenating and relaxing the feet.

Can LX-Slippers be worn outside?

No, LX-Slippers are designed for indoor use only. They are meant to provide a relaxing experience at home.

How long should I wear LX-Slippers for each session?

It is recommended to wear LX-Slippers for 10-15 minutes per session. You can customize the therapy duration according to your preference.

Is there a warranty for LX-Slippers?

Yes, LX-Slippers come with a 1-year warranty. We stand by the quality of our product and offer this warranty for your satisfaction and peace of mind.

Are LX-Slippers suitable for all foot sizes?

Yes, LX-Slippers are designed to be adjustable and fit most foot sizes comfortably.

Can LX-Slippers be used by individuals with specific foot conditions?

LX-Slippers are generally safe for use, but if you have specific foot conditions or concerns, it's advisable to consult a healthcare professional before using them.

10. Warranty:

For additional information, feedback, and warranty registration, please don't hesitate to contact us. Our dedicated support team is here to assist you.

Contact Information

Email: Contact@awastore.us

Warranty Details:

Enjoy peace of mind with a 1-year warranty for your LX-Slippers.

Warranty Registration:

To register your product for warranty coverage and access exclusive benefits, simply scan the QR code below or visit our website.



We value your satisfaction and are committed to providing exceptional customer support. Your feedback and warranty registration help us ensure your product's continued performance.

11. Disclaimer:

This device has not been evaluated by the U.S. Food & Drug Administration (FDA). Our products are designed to maintain or encourage a healthy lifestyle and are not intended to diagnose, treat, cure, or prevent any disease. Our products are low-risk, general wellness/fitness products that do not require FDA clearance, in accordance with the "General Wellness: Policy on Low-Risk Devices" draft released on January 20, 2015.

12. Contact Us:

You can contact the AWA team via Email: contact@awastore.us

You can visit www.awastore.us for more information.



American Wellness Authority



American Wellness Authority



[awa_products](#)



[awastore.us](#)



[awastore_us](#)



[awastore_us](#)